

the new



CONTOURA

by INMODE



شفا تجهیز

WWW.SHAFATAJHIZ.IR

INNOVATE
with INMODE



BODY *fx* | MINI *fx*
PLUS

Full Body Contouring Experience

GET THE ULTIMATE IN RESHAPING WITH A POWERFUL CONTOURING TOOL

BODY *fx*

The first and only non-invasive body contouring technology proven in peer reviewed and published human studies to both permanently kill adipose tissue and contract the skin.

MINI *fx*

Perfectly sized to address problematic fatty tissue in smaller areas like the arm, neck and knees.

PLUS

The first auto-adjusting thermal skin treatment for collagen remodeling and skin contraction all areas of the body.

Specifications

PLUS



RF Output Frequency:

Up to 65 W

Output Frequency:


1 MHz

Temperature Cut-Off:


Real time temperature measurements set by operator up to 43°C




DR. B. RIDENOUR

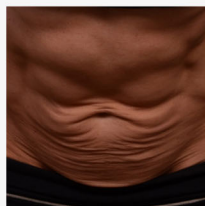

DR. B. RIDENOUR

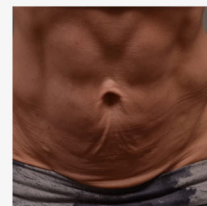

DR. S. MULHOLLAND

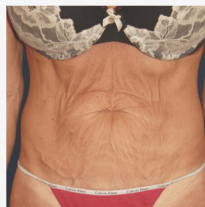

F. ORMONDE, MA




Specifications



BODYfx





Chamber Size:	2 inches x 1 inch
RF Output Frequency:	Up to 50 W
Output Frequency:	1 MHz
Negative Pressure:	Up to 500 mBar
Temperature Cut-Off:	Real time temperature monitoring and cut-off. Operator selectable up to 43°C.

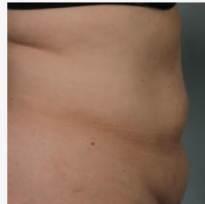



DR. D. DUNCAN

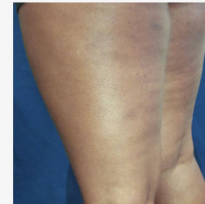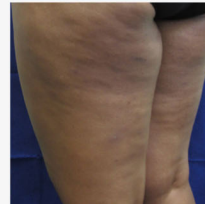

DR. J. DIAMOND

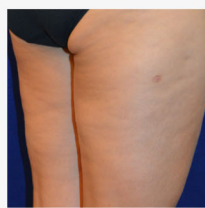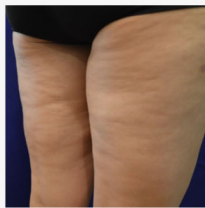

DR. J. DIAMOND

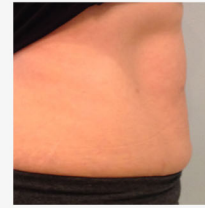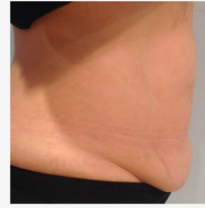


SERENITY CREEK MED SPA
BodyFX & Forma

Specifications

MINIfx



Chamber Size:	1 inch x 1 inch
RF Output Frequency:	Up to 50 W
Output Frequency:	1 MHz
Negative Pressure:	Up to 500 mBar
Temperature Cut-Off:	Real time temperature monitoring and cut-off. Operator selectable up to 43°C.

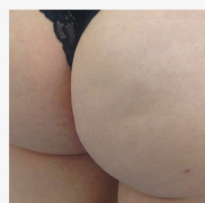


INMODE



DR. J. HELLMAN

LATEST UPDATE OF BODYFX ON APOPTOSIS

in LASERS IN SURGERY & MEDICINE JOURNAL

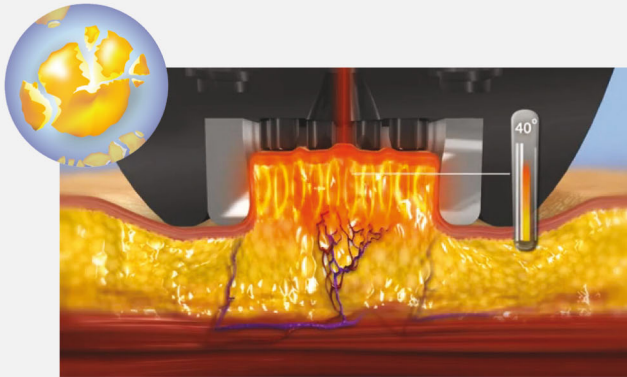

- » Fat cell apoptosis: 30% of adipocytes, in a 2.5 cm treatment zone, stained positive for APAF-1, a validated marker for apoptosis
- » Significant increase in collagen synthesis (neocollagenesis) representing an average increase of 13.7%

BENEFITS TO YOUR PRACTICE

- » Contoura has zero consumable costs or pulse counters. Unlike other body contouring devices (which may carry high per procedure cost), Contoura is a business solution that allows you to grow your patient base without having to worry about ongoing fees.
- » The use of radio frequency allows owners to provide both skin contraction and fat treatments, of all body areas, in one device. This provides multiple revenue streams to a practice's bottom line.



TOLL FREE {855} 411 2639 | WWW.INMODEMD.COM



CONTOURA IS THE RESHAPING SOLUTION TO MAXIMIZE PATIENT RESULTS

Body contouring is one of the most popular aesthetic treatments. In many cases, patients are looking for improvements in body shape without surgical intervention. The Contoura is a non-invasive body contouring solution that will provide the ultimate in lunchtime procedures.

Contoura addresses:

- » Body Reshaping
- » Cellulite Reduction
- » Skin Contraction

The Contoura is the ONLY non-invasive body contouring device proven in peer reviewed and published human studies to both permanently kill adipose tissue and contract the skin.

Provide your patients with clear and clinically proven results, so they can have the beautiful curves they want, from top to bottom.



شفا تجهیزات

WWW.SHAFATAJHIZ.IR

A.C.E. Technology

*Acquire. Control. Extend.
Technology to safely
optimize peak energy*

» ACQUIRE:

Smart RF thermistors built into the handpiece reads skin temperature 1000x per second to “ACQUIRE” therapeutic skin contraction temperatures without overheating. Measurements are clearly displayed on the system’s control screen and audible indicators alert the operator to critical changes in the temperature and impedance profile in real time.

» CONTROL:

Unlike ANY other thermal skin enhancement technology to date, the A.C.E. “CONTROL” feature is a one-of-a-kind software algorithm that allows the clinician unprecedented safety and control by providing programmable “cut-off” temperature between 35°-42°C.

» EXTEND:

Clinical evidence indicates that prolonged exposure at (41°C - 42°C) is advantageous for skin tightening. With the A.C.E. “EXTEND” feature you can now expose your patients skin, adipose tissue, or both to clinically therapeutic temperatures for as long as you want without the fear of adverse incidents in order to maximize your patients’ outcomes and drastically reduce or eliminate non-responders.



شفا تجهیز

WWW.SHAFATAJHIZ.IR